**LCD-MIKROSKOP TOUCH**  
**LCD MICROSCOPE TOUCH**

Art. No. 52-01002

- DE BEDIENUNGSANLEITUNG
- GB INSTRUCTION MANUAL
- FR MODE D'EMPLOI
- ES INSTRUCCIONES DE USO

Fig. 1

Fig. 2





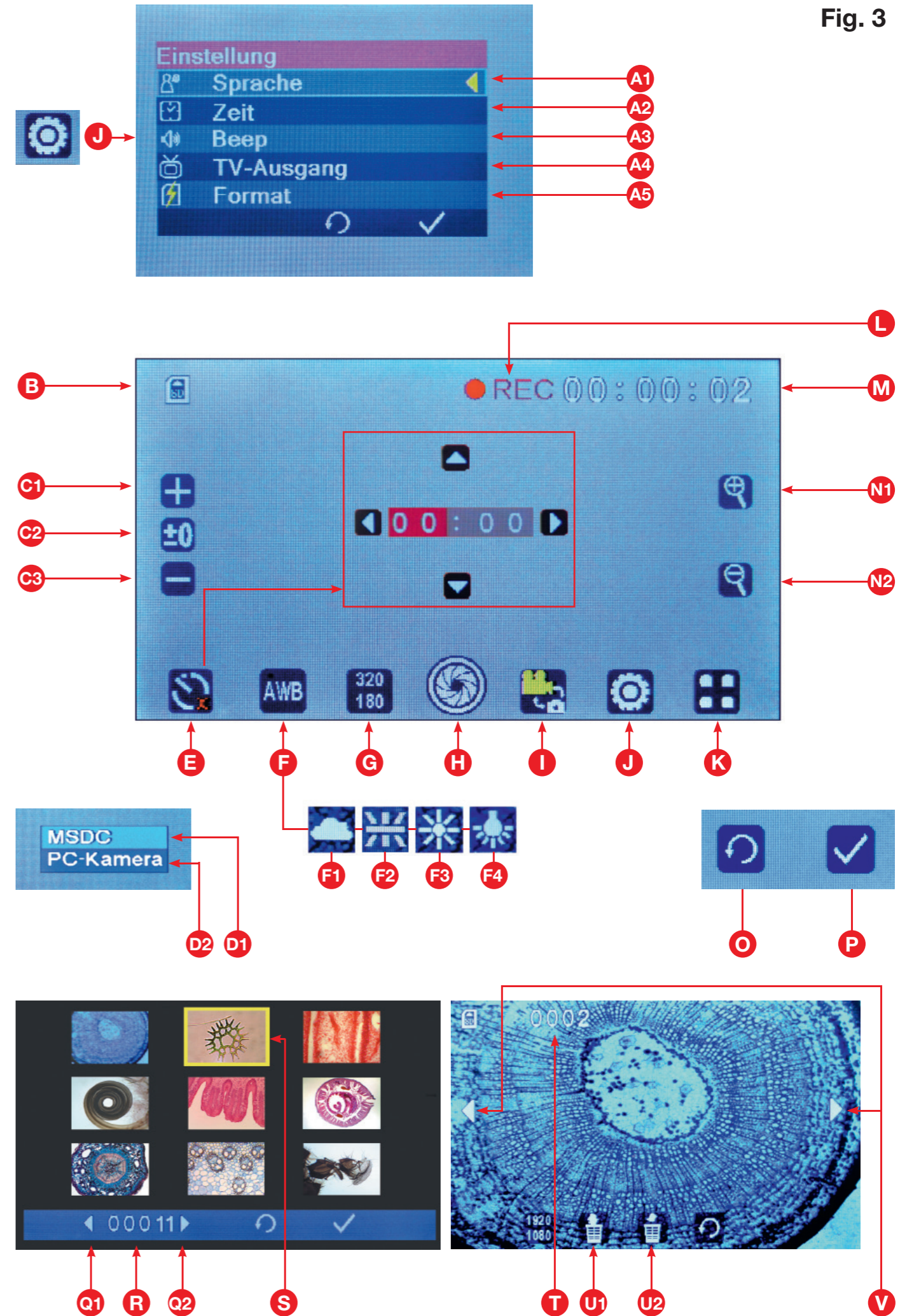
Fig. 3

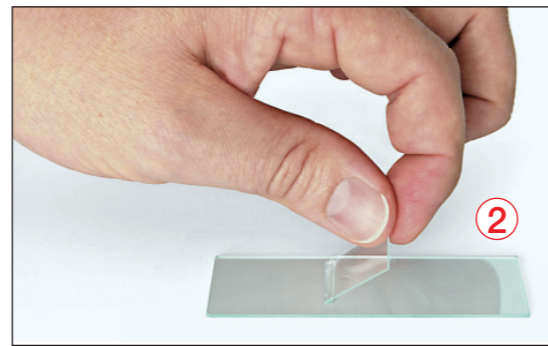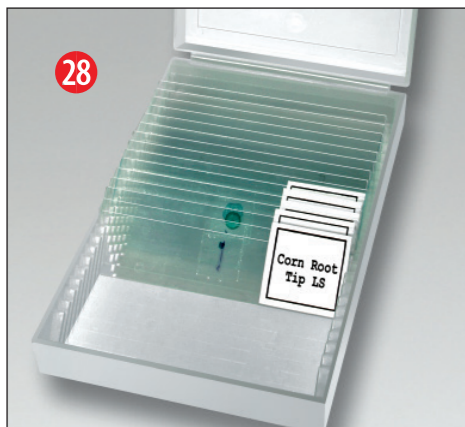**B**

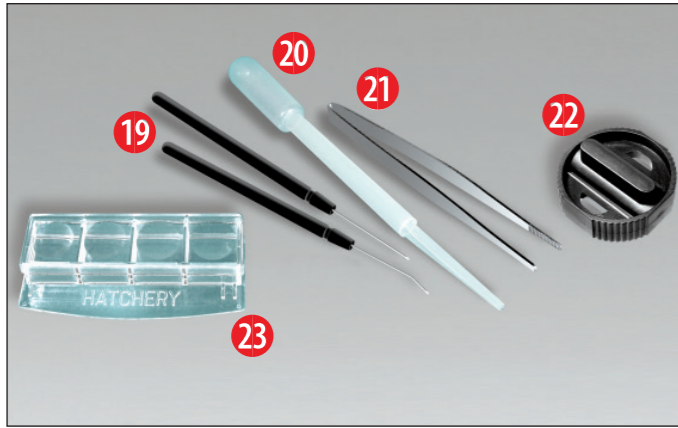optional  
optionnel

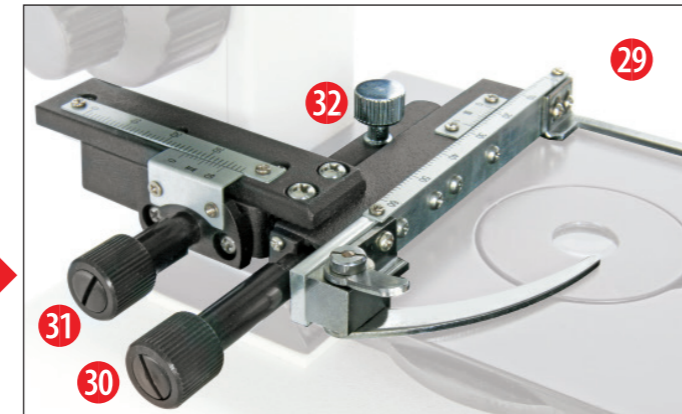DE	BEDIENUNGSANLEITUNG.....	5
GB	INSTRUCTION MANUAL.....	10
FR	MODE D'EMPLOI .....	14
ES	INSTRUCCIONES DE USO.....	18

**GARANTIE & SERVICE / WARRANTY & SERVICE**  
**GARANTIE ET SERVICE / GARANTÍA Y SERVICIO**  
**GARANZIA E ASSISTENZA / ГАРАНТИЯ.....22**

## Allgemeine Informationen

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## General Information

### About this Instruction Manual

Please read the safety instructions in this manual carefully. To avoid damage to the unit and the risk of injury, only use this product as described in the manual. Keep the instruction manual handy so that you can easily look up information on all the functions.

#### DANGER!

You will find this symbol in front of every section of text which deals with the risk of severe injury or even death in the event of improper use.

#### CAREFUL!

You will find this symbol before every section of text that deals with the risk of minor to severe injuries resulting from improper use.

#### NOTE!

You will find this symbol in front of every section of text which deals with the risk of damage to property or the environment.

### Intended Use

This product is intended only for private use. It was developed for the magnified display of things in nature.

### General Warning

#### DANGER!

Tools with sharp edges and points are often used when working with this device. For this reason, store this device and all accessories and tools in a location that is out of the reach of children. There is a RISK OF INJURY!

#### DANGER!

This device contains electronic components which operate via a power source (power supply and/or batteries). Do not leave children unattended when using the device. Only use the device as described in the manual, otherwise you run the RISK of an ELECTRIC SHOCK.

#### DANGER!

Do not expose the device to high temperatures. Use only the supplied power adapter. Do not short circuit the device or batteries, or throw them into a fire. Excessive heat or improper handling could trigger a short circuit, fires or even explosions.

#### DANGER!

Never bend, pinch, pull or run over the power and connecting cables or extensions/adapters. Protect the cables from sharp edges and heat. Before operation, check the device, cables and connections for damage.

Never use a damaged unit or a unit with damaged power cables, etc. Damaged parts must be exchanged by an authorized service centre immediately.

#### DANGER!

Children should only use the device under adult supervision. Keep packaging material (plastic bags, rubber bands, etc.) out of the reach of children! There is a RISK OF CHOKING!

#### CAREFUL!

Children must not have access to the included chemicals and liquids. Do not drink the chemicals. Wash hands thoroughly with running water after use. In the event of

contact with the eyes or mouth, rinse thoroughly with water. In the event of pain, contact a doctor immediately and take the substances with you.

#### NOTE!

Do not disassemble the device. In the event of a defect, please contact the Service Centre in your country (see the warranty card).

Do not expose the device to temperatures above 45°C.

Fig. 3

### All Touch Control Elements of the Integrated Software

The microscope is equipped with a pressure-sensitive monitor (touchscreen). A gentle press on the selected control element is sufficient and the corresponding sub-menu for the desired setting is displayed. The left side shows all the control menus of the integrated microscope software.

- (A1) Language selection
- (A2) Time setting, time on/off
- (A3) Audio signal on/off
- (A4) NTSC/PAL TV output selection
- (A5) Format SD card (max: 32 GB)
  
- (B) SD card in device,  
No SD card in device, symbol crossed out in red
  
- (C1) Increase the sensitivity of the installed camera
- (C2) Indication of the changed sensitivity value
- (C3) Reduce the sensitivity of the installed camera
  
- (D1) If connected to a PC via USB cable:  
data transfer to the computer activated
- (D2) If connected to a PC via USB cable:  
image transmission to the computer activated
  
- (E) Video recording time setting:  
left/right arrows: minutes/seconds selection  
Up/down arrows: minutes/seconds setting
  
- (F) Adjust the image colour temperature (calibrate image)  
AWB setting = automatic calibration  
Change the automatic calibration values:
  
- (G) Set the resolution for photos or video  

Photos:	Video:
640 x 360 pixels	320 x 180 pixels
1920 x 1080 pixels	640 x 360 pixels
2048 x 1152 pixels	
2560 x 1440 pixels	
3072 x 1728 pixels	
3648 x 2048 pixels	

- (H) Capture a photo or start/stop video recording
- (I) Switch between photo capture and video recording
- (J) Basic program settings
- (K) Complete overview of captured images
- (L) Video recording active
- (M) Length of video recording
  
- (N1) Enlarge image
- (N2) Reduce image
  
- (O) Return to the overview
- (P) Confirm the selected setting

(Q1/Q2) Select image number

- (R) Indication of the selected image number
- (S) Yellow frame marks the selected image
- (T) No. of the image currently displayed on the monitor

- (U1) Delete all images
- (U2) Delete a single image marked with a yellow frame

(V) Arrows for selecting the previous or next image

#### Note:

If no control elements are displayed, a light press on the screen is sufficient to make all the basic control elements reappear!

### All parts (Fig. 1)

- 1 Monitor
- 2 Monitor on/off button
- 3 Revolving turret
- 4 Interchangeable objectives with varying magnification
- 5 Microscope table
- 6 LED transmitted light illumination
- 7 Power supply
- 8 Colour filter wheel with pinhole apertures
- 9 Storage compartment for small items
- 10 Focusing wheel
- 11 Rear power connection

### All parts (Fig. 2)

- 12 SD card slot, max. 32 GB
- 13 Mini USB connection
- 14 AV output
- 15 LED incident light illumination
- 16 Microscope slide clips
- 17 On/off switch with integrated dimmer for incident light illumination
- 18 On/off switch with integrated dimmer for transmitted light illumination

### All parts (Fig. 3)

- 19 Dissecting needles
- 20 Pipette
- 21 Tweezers
- 22 MicroCut forceps
- 23 Shrimp hatchery
- 24 Yeast
- 25 Gum media
- 26 Sea salt
- 27 Shrimp eggs
- 28 Box containing microscope slides, cover slips and permanent specimens
- 29 Cross table (optional)
- 30 Transverse adjusting screw
- 31 Longitudinal adjusting screw
- 32 Fastening screw

### Scope of delivery:

- LCD touch microscope with slide clips and colour filter wheel with pinhole apertures
- Two dissecting needle
- Pipette
- Tweezers
- MicroCut forceps
- Yeast, gum media, sea salt, shrimp eggs, shrimp hatchery
- Box containing microscope slides, cover slips and permanent specimens
- Mini USB connection cable
- Power supply
- Software CD

### General/location

Select a suitable location before starting to set up your microscope. First make sure that your microscope is standing on a stable, vibration-free surface. In order to put the microscope into operation, the power supply provided must be connected to the microscope (Fig. 1, image 7/11) and a 230 V plug socket. The appropriate illumination is switched on using the switch (Fig. 1, image 17/18), and the brightness is controlled for optimum illumination of the specimen by turning the selection wheel. (For safety reasons, please only use the power supply provided. Other power supplies may not meet the required technical specifications. Damage to the device caused by other power supplies is not covered by the guarantee.)

## 1. Installation and use of the software

### 1.1. Software and installation information

A software CD is included with your microscope. The software and drivers on this CD must be installed on your computer in order to use your MicOcular. Once you connect the MicOcular to your computer after installation, you can view pictures on and save them to your computer using Photomizer. To install the software and driver correctly, simply follow the installation steps.

### 1.2. Software installation

1. Important: Before inserting the CD, first plug the USB cable into the USB port on your PC. Windows will now recognize the new device, and will indicate this in a notification window. Now please click on "Cancel", since we will have the driver for the device automatically installed later.
2. Now insert the CD-ROM that came in your package into the CD/DVD drive of your computer. The installation menu starts automatically. If it does not, go to the Windows explorer and select the CD/DVD drive (most of the time, it's the "D" drive, but it can have another letter). From there, start the file "menu.exe" by double-clicking with the left mouse button.

### 1.2.1. Driver installation

#### CAUTION!

Save all current data and close all open programs.

To install the driver software, click on the menu point "Install Driver" with the mouse cursor. Then follow the installation program instructions. During the software installation the correct driver for your operating system will be installed automatically. No manual input is needed. In rare cases the device may not be recognised by your computer. As a rule you need then only install the driver manually from the CD. If this fails please refer to the troubleshooting chapter that follows.

Restart the computer after installing the driver.

### 1.2.2. Installation of the Image Editing Software Photomizer

The image editing software "Photomizer" is located on the software CD. You can edit your pictures here.

1. To install the software, click once on the menu point "Install Photomizer" with the left mouse button.
2. The Photomizer Software requires Microsoft .NET Framework 3.5, which it will install if it is not already on your system. If it is already on your system, you may skip to step 5.
3. In the welcome window, please accept the Microsoft license agreement, and then click "Install". The installation can take a few minutes.
4. As soon as everything is installed, click "Finish".
5. Now you will be presented with a choice, in which you can choose your language.  
Make your selection and confirm it by clicking „OK“.
6. When you see "Welcome", click on „Next“.
7. In the next window, you will be asked for the „Destination Folder“. Here, just click on „Next“.
8. The „Ready to install Photomizer“ window appears. Select whether a shortcut should be created on the desktop or the quick start bar.
9. Now, the window with the Setup-Status will appear – here, a progress bar will inform you about the ongoing installation. This process can take a few minutes.
10. The window "Photomizer is being installed" appears. Click on "Finish". The installation ends.

### 1.3. Troubleshooting: the device isn't recognised by the computer (USB)

Please contact our Service: [service@bresser.de](mailto:service@bresser.de)

### 2. Electric LED illumination with dimmer

Transmitted light illumination (Fig. 1, image 6) is used for transparent specimens (specimens on microscope slides). Select incident light illumination (Fig. 2, image 15) to observe opaque objects. The simultaneous use of transmitted and incident light illumination is only advisable with semi-transparent objects. This mode is not recommended for transmitted light objects on slides since too much light may be reflected on the slide.

### 3. Colour filter/aperture combi wheel

The colour filter/aperture combi wheel (Fig. 1, image 8) beneath the microscope table can be used to influence the image quality when observing transparent specimens. Fine details are displayed better depending on the colour and object. A coloured lower light in combination with a white upper light can also optimise image details when observing objects under incident light (e.g. transparent objects).

Depending on the aperture used, the corresponding light bundling can also influence the sharpness, depth of field, contrast and detail resolution during observation.

### 4. Observation

Turn the nosepiece (Fig. 1, image 3) to select the 4x objective for observing objects. Depending on the object, select upper or lower lighting for optimum illumination. Place a specimen on the microscope table (Fig. 1, image 5) directly under the objective and secure it using the slide clips (Fig. 2, image 16).

**Important:** The object to be observed must lie precisely in the beam path of the transmitted light illumination. Now slowly turn the focusing wheel (Fig. 2, image 10) until the image on the monitor is reproduced sharply. The dimmers (Fig. 2, image 17) and (Fig. 2, image 18) can be used to optimally adjust the contrast and quality of the image of the object under observation.

Turn the nosepiece (Fig. 1, image 3) to switch objectives if increased magnification is desired. Incidentally: The higher the level of magnification, the more light is required for good image quality.

Incident light illumination is optimised for the 4x objective. The specimen is generally no longer fully illuminated if the 10x or 40x objective is used.

### 5. Viewed Object – condition and preparation

Both transparent and non-transparent specimens can be examined with this microscope, which is a direct as well as transmitted light model. If opaque specimens are examined – such as small animals, plant parts, tissue, stone and so on – the light is reflected from the specimen through the lens and eyepiece, where it is magnified, to the eye (reflected light principle, switch position I). If opaque specimens are examined the light from below goes through the specimen, lens and eyepiece to the eye and is magnified en route (direct light principle, switch position II). Many small organisms of the water, plant parts and finest animal components have now from nature these transparent characteristic, other ones must be accordingly prepared. Is it that we make it by means of a pre-treatment or penetration with suitable materials (media) transparent or thus that we cut finest wafers off of them (hand cut, MicroCut) and these then examine. With these methods will us the following part make familiar.

### 6. Creation of thin preparation cuts

Specimens should be sliced as thin as possible, as stated before. A little wax or paraffin is needed to achieve the best results. A candle can be used for the purpose. The wax is put in a bowl and heated over a flame.

#### DANGER!

Be extremely careful when working with hot wax. There is a risk of burns.

The specimen is then dipped several times in the liquid wax. The wax is finally allowed to harden. Use a MicroCut or knife/scalpel (carefully) to make very thin slices of the object in its wax casing.

#### DANGER!

Be extremely careful when using the knives/scalpels or the MicroCut. There is an increased risk of injury due to the sharp edges!

These slices are then laid on a glass slide and covered with another.

### 7. Creation of an own preparation

Put the object which shall be observed on a glass slide and give with a pipette (Fig. B 1) a drop of distilled water on the object (Fig. B 2).

Set a cover glass (in each well sorted hobby shop available) perpendicularly at the edge of the water drop, so that the water runs along the cover glass edge. Lower now the cover glass slowly over the water drop.

### 8. Gum Media

Gum media can be used to embed or preserve specimens. Warm the gum media slightly in a water bed to make it easier to use. It is used similarly to how water is used when preparing a specimen (image B1), except that a dissecting needle is used in this case. Vaseline or a clear nail varnish, for example, can also be used to preserve specimens.

### 9. Connection to a PC/laptop

Connect the USB cable to the socket (Fig. 1, image 13) on the microscope and a free USB port on your PC.

Press „MSCD“ to transfer images directly from the SD card.

Press „PC camera“ to transmit the live image from the camera. The resolution corresponds to the video resolution (Fig. 3,

image G). In order to transmit a live image, it is necessary to install the drivers and software contained on the CD provided.

### 10. The lifecycle of the saltwater prawn

The saltwater prawn or artimia salina to scientists has an unusual and interesting lifecycle. The female's eggs are hatched without any male shrimp having to fertilise them. The resultant baby prawns are all female. Under unusual circumstances such as when a swamp is drained the eggs may produce male prawns. These males fertilise the female's eggs, resulting in a specific type of eggs. These are called winter eggs and have a thick shell as protection. They're pretty rugged and can survive the swamp or lake drying out causing the death of the entire prawn population for up to a decade in a form of hibernation. The eggs hatch once the right ambient conditions again obtain. The eggs supplied (Fig B 27) are of this type.

### 11. Hatching of the salt water prawn

To hatch the prawn it is essential to first have a saline solution suited to the prawn's needs. Fill half a litre of rain- or fresh-water in a container. Let it stand for about thirty hours. As water evaporates over time it's a good idea to have a second container of such water left standing for thirty-six hours. Once it's stood for this length of time pour half of the sea salt supplied into one of the containers (Fig B 26) and stir until it has dissolved. Then pour some of it into the prawn breeding plant (Fig B 23). Add a few eggs and close the lid. Put it somewhere with plenty of light but not in the direct sun. The temperature should be approximately 25° C. The prawns will hatch in two or three days at this temperature. Should any water evaporate during this time replace it from the second container.

### 12. The saltwater prawn under the microscope

What comes out of the egg is known as a nauplius larva. Use the pipette (Fig B 20) to put some of them on a slide for examination. They will move in the solution using their hair like limbs. Remove a few daily from the container for examination under the microscope. If you do so and save the pictures made with the MicOcular you will then have a seamless record of the prawn's lifecycle. You can remove the upper lid of the prawn bleeding plant and put the whole thing under the microscope. The larvae will mature in six to ten weeks depending on ambient temperature. You will soon have bred an entire generation of saltwater prawns that constantly reproduce.

### 13. Notes on cleaning

#### 1. Avoid disassembling the microscope!

Because of being a precise optical instrument, do not disassemble the microscope. That may cause serious damage to its performance. Do not disassemble the objectives.

#### 2. Avoid contamination!

You should make sure that dust or moisture do not come in contact with your microscope. Avoid putting fingerprints on any optical surfaces! If dust or dirt nevertheless have gotten onto your microscope, or the accessories, remove it as described below.

#### 3. Cleaning non-optical parts:

Disconnect the plug from mains socket before cleaning non-optical parts of the microscope (e. g. the frame or the stage). Firstly, remove dust or dirt from your microscope, or accessories, using a soft brush. Then clean the soiled area with a soft, lint-free cloth; the cloth may be moistened with a mild detergent solution. Be sure that the microscope is dry before using. For cleaning the painted parts or plastic components, do NOT use organic solvents such as alcohol, ether, acetone, xylene or other thinner etc.!

### 4. Cleaning optical parts:

In order to enhance the optical quality, the lenses of eyepieces or objectives are coated. They should not be wiped because dry dirt or dust may scratch the coating. It is best to remove the parts to be cleaned from the frame prior to cleaning. Always blow loose dust away first. Use lens tissue of good quality or soft fabric moistened with a small amount of pure alcohol (available in the chemist's shop); wipe the lens surface clean.

### 5. When not in use:

After usage, cover the microscope with the dust cap (made of plastic) and place it where it is dry and mouldless. In the case of non-use for a longer time, you should stow the microscope and the accessories in their correct containers again. We suggest the storage of all objectives and eyepieces in a closed container with drying agent.

**Remember: A well maintained microscope will keep its optical quality for years and thus maintain its value.**

### DISPOSAL

Dispose of the packaging materials properly, according to their type (paper, cardboard, etc). Contact your local waste disposal service or environmental authority for information on the proper disposal.

Do not dispose of electronic devices in the household garbage! As per the Directive 2002/96/EC of the European Parliament on waste electrical and electronic equipment and its adaptation into national law, used electronic devices must be collected separately and recycled in an environmentally friendly manner.

Please take the current legal regulations into account when disposing of your device. You can get more information on the proper disposal from your local waste disposal service or environmental authority.

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#### **GB WARRANTY & SERVICE**

The regular guarantee period is 2 years and begins on the day of purchase. To benefit from an extended voluntary guarantee period as stated on the gift box, registration on our website is required. You can consult the full guarantee terms as well as information on extending the guarantee period and details of our services at [www.bresser.de/warranty\\_terms](http://www.bresser.de/warranty_terms). Would you like detailed instructions for this product in a particular language? Then visit our website via the link below (QR code) for available versions. Alternatively you can also send an email to [manuals@bresser.de](mailto:manuals@bresser.de) or leave a message on +49 (0) 28 72 - 80 74-220\*. Please always state your name, precise address, a valid phone number and email address, as well as the article number and name.

\*Number charged at local rates in Germany (the amount you will be charged per phone call will depend on the tariff of your phone provider); calls from abroad will involve higher costs.

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**Bresser GmbH**